

How to Check a Blood Pressure Manually

Purpose

This standard operating procedure describes how to accurately measure a patient's blood pressure using a manual cuff (sphygmomanometer) and stethoscope. It covers patient preparation, cuff placement, estimation of systolic pressure to avoid the auscultatory gap, auscultation and recording of readings, and interpretation of results according to current clinical guidelines.

Scope

This procedure applies to any clinical staff member responsible for performing manual blood pressure measurements on patients in a seated or bedside setting.

Required equipment

- Stethoscope
- Manual blood pressure cuff (sphygmomanometer), correctly sized for the patient's arm
- Hand sanitizer or handwashing facilities
- Documentation materials to record the reading

Procedure

Step 1: Perform hand hygiene

Wash your hands thoroughly or use hand sanitizer before beginning the procedure.

Step 2: Gather your supplies

Obtain a stethoscope and a manual blood pressure cuff before approaching the patient.

Step 3: Prepare the patient

- Ensure the patient is sitting down.
- Position the patient's arm at heart level, resting comfortably on a table or armrest.
- Make sure the patient's legs are uncrossed.

- If the patient is lying in bed, adjust the arm so it remains at heart level.



Nurse and patient seated at a table, with the patient's arm at heart level while the nurse explains positioning

Step 4: Select the correct cuff size

Choose a blood pressure cuff that fits the patient's arm properly. Safety note: Avoid using a cuff that is too large or too small, as this can result in inaccurate readings.



Nurse explaining the importance of correct cuff size, gesturing with both hands, with the cuff visible on the table

Step 5: Identify the brachial artery

Palpate (feel for) the brachial artery, located in the bend of the arm (antecubital fossa). This is the artery you will use to measure blood pressure.



Nurse demonstrating the location of the brachial artery on the patient's exposed arm, pointing to the upper arm

Step 6: Understand the blood pressure readings

- The systolic number (top number) is the first sound you hear as you release air from the cuff.

- The diastolic number (bottom number) is the point where you no longer hear the sound.




Nurse explaining the systolic number, with text overlay: "Systolic Number (Top Number/First Sound)"



Nurse explaining the diastolic number, with text overlay: "Diastolic Number (Bottom Number/Last Sound)"

Step 7: Apply the blood pressure cuff

Wrap the cuff snugly around the patient's upper arm, aligning the artery marker (arrow) on the cuff with the brachial artery. Ensure the gauge and tubing are positioned for easy access and visibility.

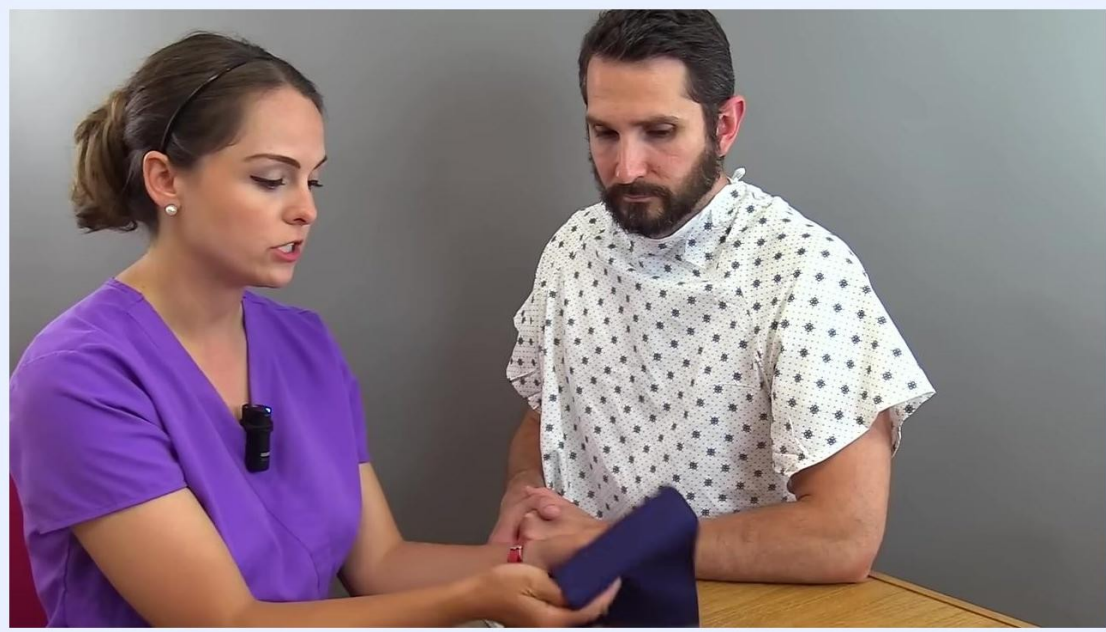

Close-up of the blood pressure cuff showing the arrows and artery marker, with the nurse pointing to the indicator

Step 8: Position the cuff correctly for measurement

- Place the blood pressure cuff on the patient's upper left arm.
- Ensure the arrow on the cuff points toward the brachial artery.

- Position the cuff approximately two inches above the bend of the arm (antecubital fossa).



Nurse positioning the cuff about two inches above the bend of the patient's left arm, with the arrow on the cuff visible and the patient holding their sleeve up

Step 9: Estimate the systolic pressure by palpation

- Palpate (feel) the brachial artery with your fingers.
- Inflate the cuff while palpating the artery.
- Continue inflating until you can no longer feel the pulse in the brachial artery.


- Observe the gauge at the exact point when the pulse disappears — this value is your estimated systolic pressure.



Nurse explaining the estimation process, with the cuff on the patient's arm and the gauge visible

Step 10: Record the estimated systolic pressure

Note the reading on the gauge at the moment the pulse disappears. For example, if the pulse disappears at 100 mmHg, record this as the estimated systolic pressure.

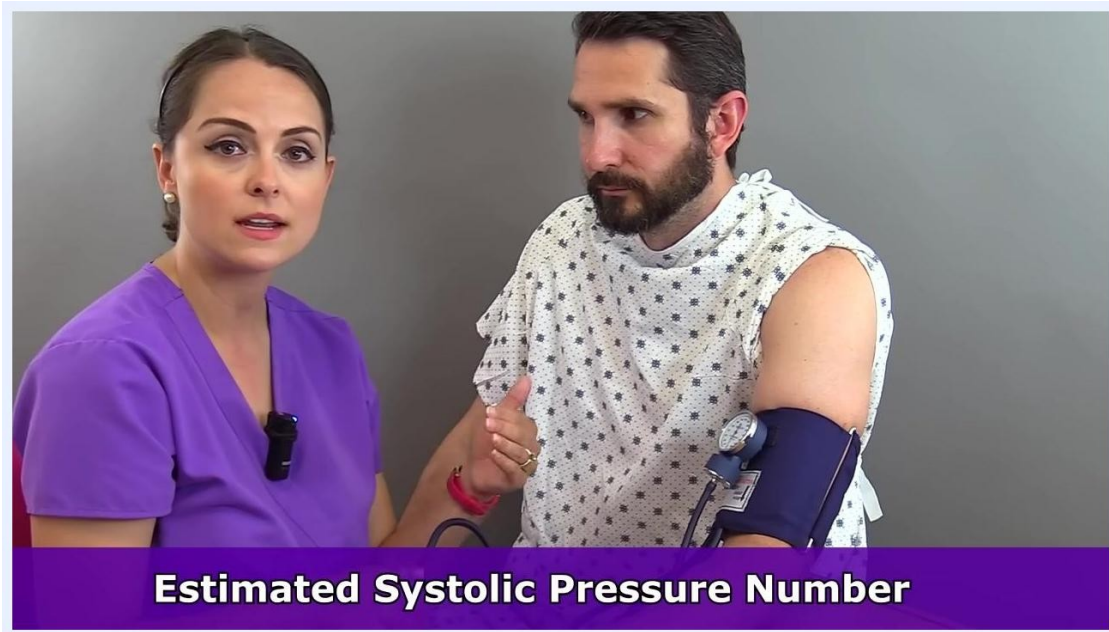

Nurse and patient with the cuff in place, nurse referencing the gauge, both focused on the process

Step 11: Deflate the cuff completely

Release all air from the cuff after noting the estimated systolic pressure. Wait 30 to 60 seconds before proceeding to the next step to allow blood flow to return to normal.


Step 12: Prepare to take the actual blood pressure reading

Re-inflate the cuff to a value 30 mmHg higher than your estimated systolic pressure. For example, if your estimated systolic pressure is 100 mmHg, inflate the cuff to 130 mmHg.

Safety note: This step helps avoid missing the auscultatory gap, a period of abnormal silence that can occur (especially in patients with hypertension) and may cause you to miss the true systolic reading.



Nurse explaining the importance of inflating 30 mmHg above the estimated systolic pressure, with text overlay: "Estimated Systolic Pressure Number"

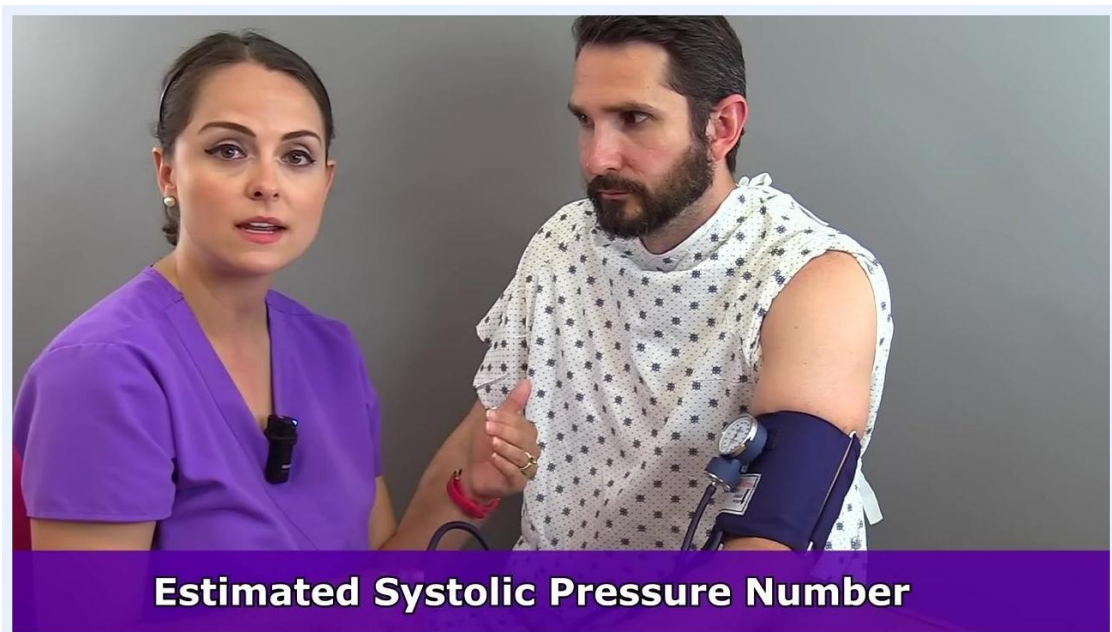

Nurse discussing the auscultatory gap, with the cuff on the patient's arm, both focused on the process

Example: If the estimated systolic pressure is 100 mmHg, inflate the cuff to 130 mmHg for the actual measurement.



Close-up of a sphygmomanometer gauge showing 100 mmHg, labeled "ReliOn", with a clear needle position

Step 13: Select the appropriate side of the stethoscope


Use either the bell or the diaphragm of your stethoscope. Prefer the bell for blood pressure measurement, as it is best at picking up low-pitched Korotkoff sounds.



Nurse holding a stethoscope, preparing to use the bell side, with the cuff on the patient's arm, both focused on the procedure

Step 14: Place the stethoscope over the brachial artery


Position the bell (or diaphragm) of the stethoscope lightly over the brachial artery, just below the cuff. Safety note: Do not press too hard, as fully compressing the artery can occlude it and prevent you from hearing the sounds.

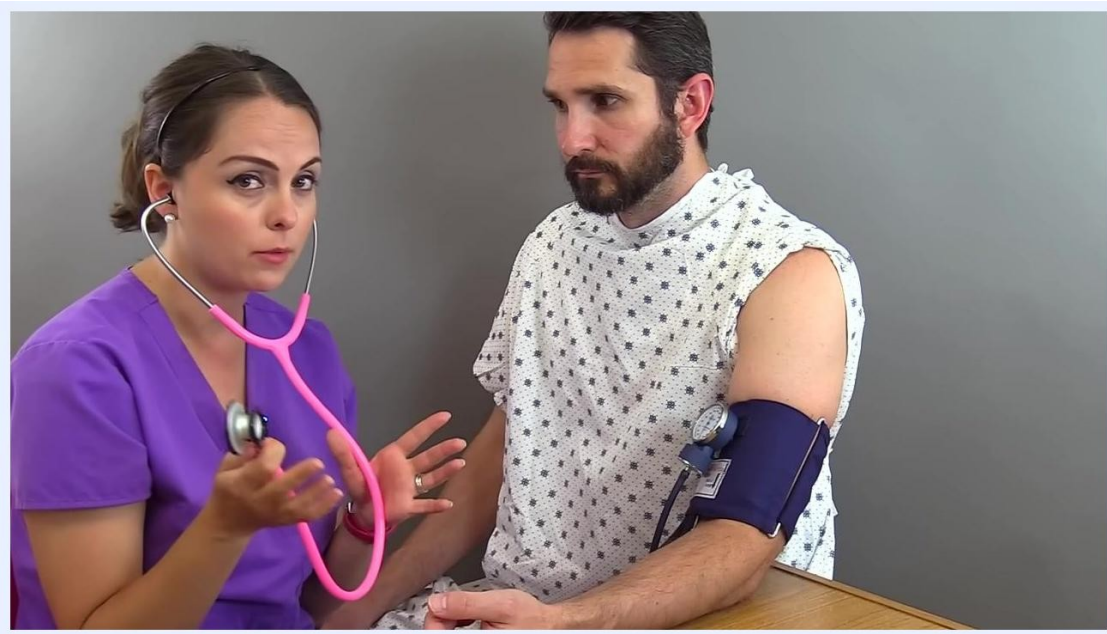

Nurse explaining stethoscope placement, holding the stethoscope near the patient's arm, cuff in place

Step 15: Inflate the cuff to the correct pressure

Inflate the cuff to 130 mmHg (30 mmHg above the previously estimated systolic pressure of 100 mmHg). Ensure the gauge reads 130 mmHg before beginning to deflate.

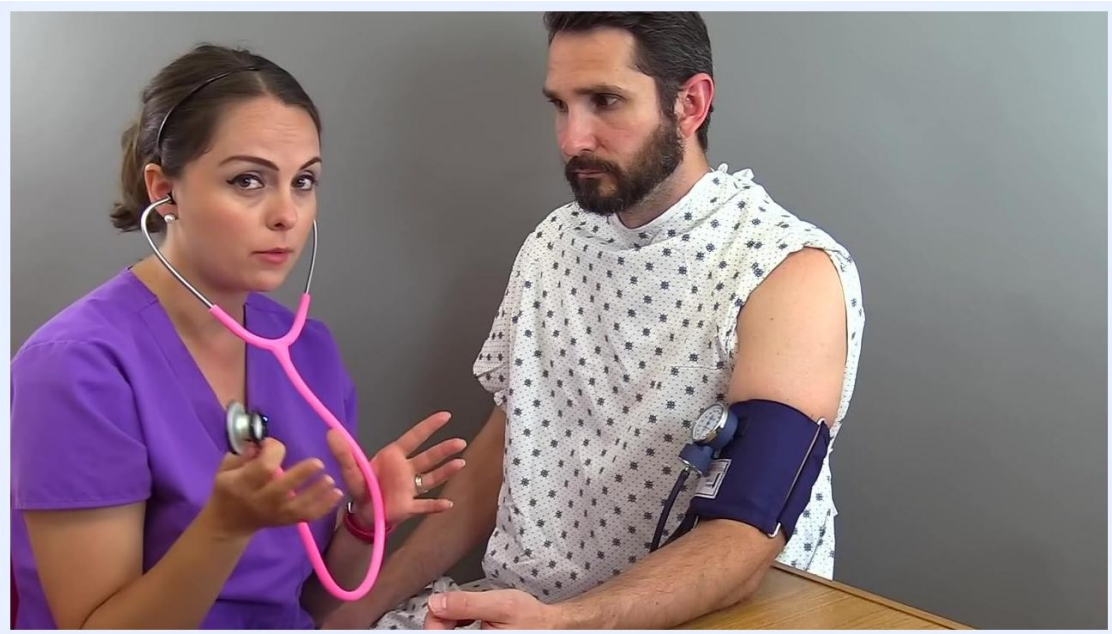

Nurse inflating the cuff, stethoscope in place, gauge visible, patient's arm positioned for measurement


Step 16: Deflate the cuff at the proper rate

Slowly release the air from the cuff at a rate of about 2 mmHg per second. Listen carefully for the Korotkoff sounds through the stethoscope.


Step 17: Identify the systolic and diastolic readings

- Note the gauge reading at the first appearance of repetitive tapping sounds (systolic pressure). In this example, the first sound is heard at 104 mmHg.

- Continue to listen as the pressure drops. Note the gauge reading at the last sound you hear (diastolic pressure). In this example, the last sound is heard at 78 mmHg.

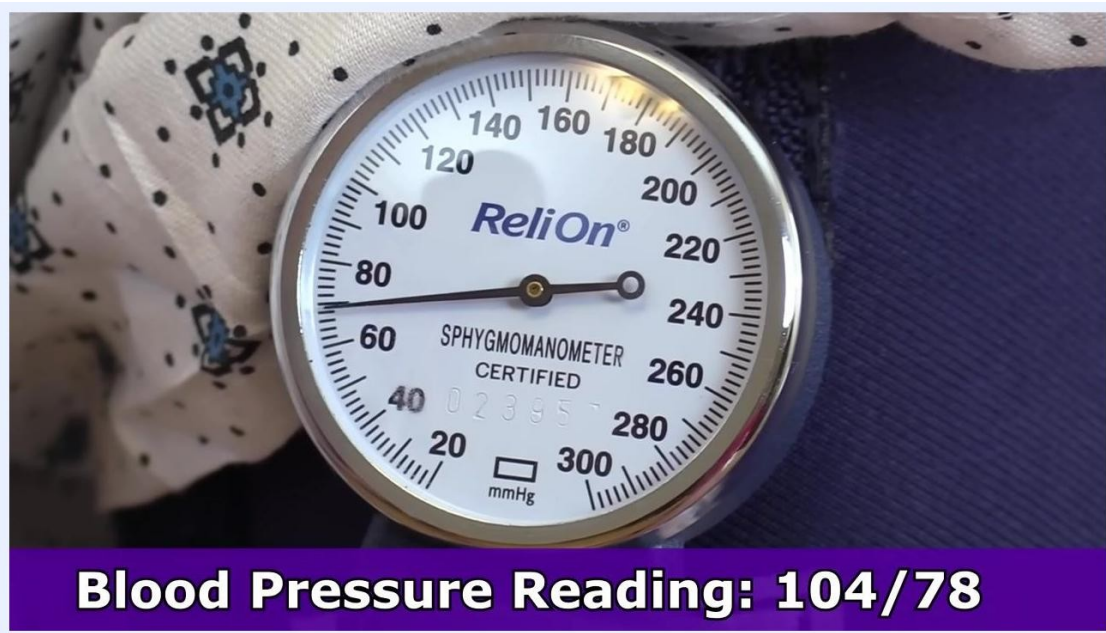


Close-up of a sphygmomanometer gauge showing the needle at 104 mmHg, labeled "ReliOn"

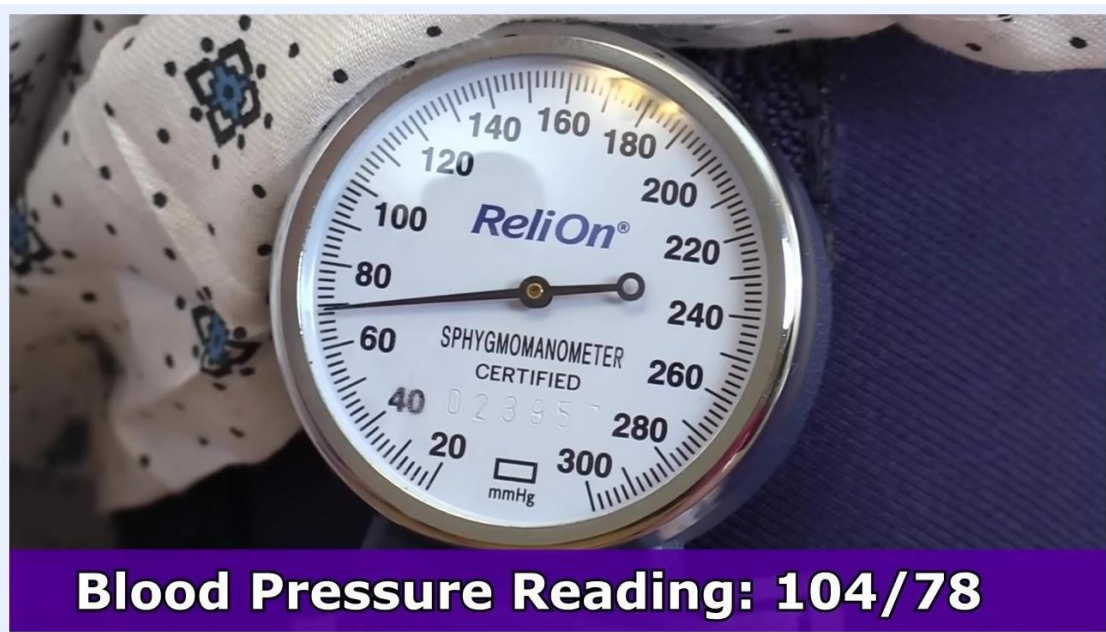

Close-up of a sphygmomanometer gauge showing the needle at 78 mmHg, labeled "ReliOn"

Step 18: Record the blood pressure reading

- Document the blood pressure as systolic over diastolic (e.g., 104/78 mmHg).
- Note which arm was used for the measurement.
- If the cuff is not disposable, clean it according to facility protocol.

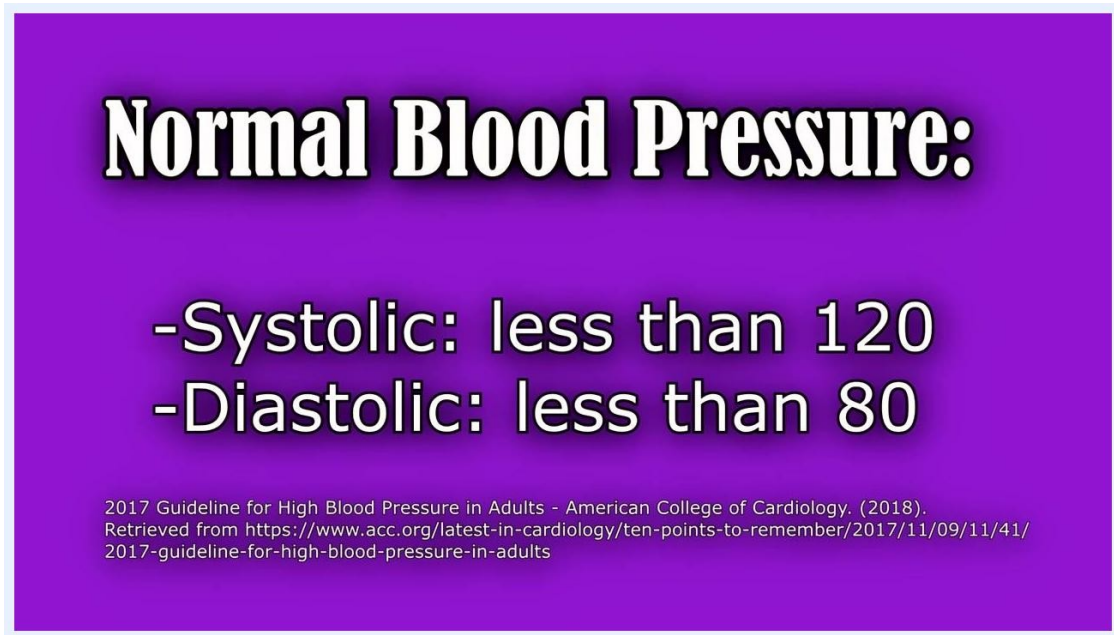

Gauge showing 78 mmHg, with text overlay: "Blood Pressure Reading: 104/78"

Step 19: Remove the cuff and stethoscope

Fully deflate the cuff and remove it from the patient's arm. Remove the stethoscope from the patient's arm.

Verification: Interpreting blood pressure categories

Use the following categories, based on the 2017 American College of Cardiology guidelines, to interpret the recorded reading:

- Normal blood pressure: Systolic less than 120 mmHg; Diastolic less than 80 mmHg.
- Elevated blood pressure: Systolic 120–129 mmHg; Diastolic less than 80 mmHg.
- Hypertension Stage 1: Systolic 130–139 mmHg; Diastolic 80–89 mmHg.

- Hypertension Stage 2: Systolic greater than 140 mmHg; Diastolic greater than 90 mmHg. See a doctor for a proper diagnosis if readings fall within this range. (Reference: 2017 Guideline for High Blood Pressure in Adults – American College of Cardiology.)

Normal Blood Pressure:

- Systolic: less than 120
- Diastolic: less than 80

2017 Guideline for High Blood Pressure in Adults - American College of Cardiology. (2018).
Retrieved from <https://www.acc.org/latest-in-cardiology/ten-points-to-remember/2017/11/09/11/41/2017-guideline-for-high-blood-pressure-in-adults>

Text slide: "Normal Blood Pressure: Systolic less than 120, Diastolic less than 80"

Elevated Blood Pressure:

- Systolic: 120–129
- Diastolic: less than 80

2017 Guideline for High Blood Pressure in Adults - American College of Cardiology. (2018).
Retrieved from <https://www.acc.org/latest-in-cardiology/ten-points-to-remember/2017/11/09/11/41/2017-guideline-for-high-blood-pressure-in-adults>

Text slide: "Elevated Blood Pressure: Systolic 120–129, Diastolic less than 80"

Hypertension Stage II:

-Systolic: Greater than 140

-Diastolic: Greater than 90

*See a doctor for a proper diagnosis

2017 Guideline for High Blood Pressure in Adults - American College of Cardiology. (2018).
Retrieved from <https://www.acc.org/latest-in-cardiology/ten-points-to-remember/2017/11/09/11/41/2017-guideline-for-high-blood-pressure-in-adults>

Text slide: "Hypertension Stage II: Systolic greater than 140, Diastolic greater than 90, See a doctor for a proper diagnosis", with ACC guideline citation

Summary

To complete the manual blood pressure check:

- Confirm you have recorded the blood pressure reading and cleaned any reusable equipment.
- Ensure the patient is comfortable and explain the results if appropriate.
- Compare the reading against the blood pressure categories above.
- Advise the patient to consult a healthcare provider for further evaluation if their reading falls into the hypertension range.

What's next

For additional educational resources, visit RegisteredNurseRN.com to access more tutorials, watch further videos, subscribe to the channel for updates on new content, and like and share videos to support the channel.

Hypertension Stage II:

- Systolic: Greater than 140
- Diastolic: Greater than 90

*See a doctor for a proper diagnosis

2017 Guideline for High Blood Pressure in Adults - American College of Cardiology. (2018).
Retrieved from <https://www.acc.org/latest-in-cardiology/ten-points-to-remember/2017/11/09/11/41/2017-guideline-for-high-blood-pressure-in-adults>

End screen with RegisteredNurseRN.com branding, a "More Videos" link, a "Subscribe" button, and a reminder to give a thumbs up and share